



tRide

Simple ... is the Best
for the Child and its Parents



“It would be great if there were a rehab stroller that combined the most important functions in one easy step!”

Frauke Bessel, Bremervörde




Three Components = Simple Solution!



Surveys show that parents want it simple.

Parents are often overwhelmed by the technical equipment of many rehab strollers. For this reason they only use some of the settings. Parents would like things to be more straightforward and the focus to be on the main functions.

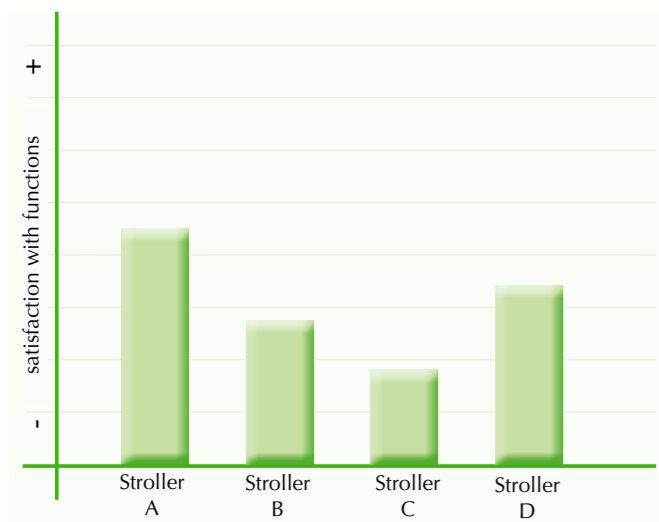

IGAP - Institut für Innovationen im Gesundheitswesen und angewandte Pflegeforschung e.V.
Source: february 2012

tRide – Easy is Better



Simply Hassle-Free



- The tRide offers lots of storage room under the seat for purchases, diapers and much more
- Practical, one-hand adjustment with dual function
 - Adjusting the angle of the seat
 - Removing the seat from the frame
- Easy to put in the trunk, saves space and is easy on your back

Relaxed Sitting – Maximum Mobility



- tRide headrest adjusts to meet the needs of your child
- Feeling of safety and comfort because of the supporting surfaces for the arms
- Soft, breathable seat cover for a cozy feeling
- Facilitates getting into and out of the buggy independently

Mobility



Sophisticated Functionality - tRide Rehab Stroller



Seat Cover
Made from durable, high quality fabric that is washable at 40° C

Footrest
Angle and height adjustable, fold up function and optional foot straps

Seat Depth and Width
Adjustable without tools using 4-star handles located below the seat even while the child is sitting in the rehab stroller

Basket or vent tray
For lots of storage, usable as required and holds up to 5 kg



Headrest
Three-dimensionally adjustable for ideal setting and support

Central Release
Using one hand the seat angle can be adjusted or the seat can be removed from the frame

Backrest Angle
Simple and easy to adjust with the central locking lever

Removing the Seat from the Frame
At a glance and in one simple step




Learn more!

All the settings and optional accessories are at www.thomashilfen.com

"A rehab stroller has to be easy and hassle-free for the parents to handle. This also ensures a good sitting position even at home."

IGAP - Institut für Innovationen im Gesundheitswesen und angewandte Pflegeforschung e.V.

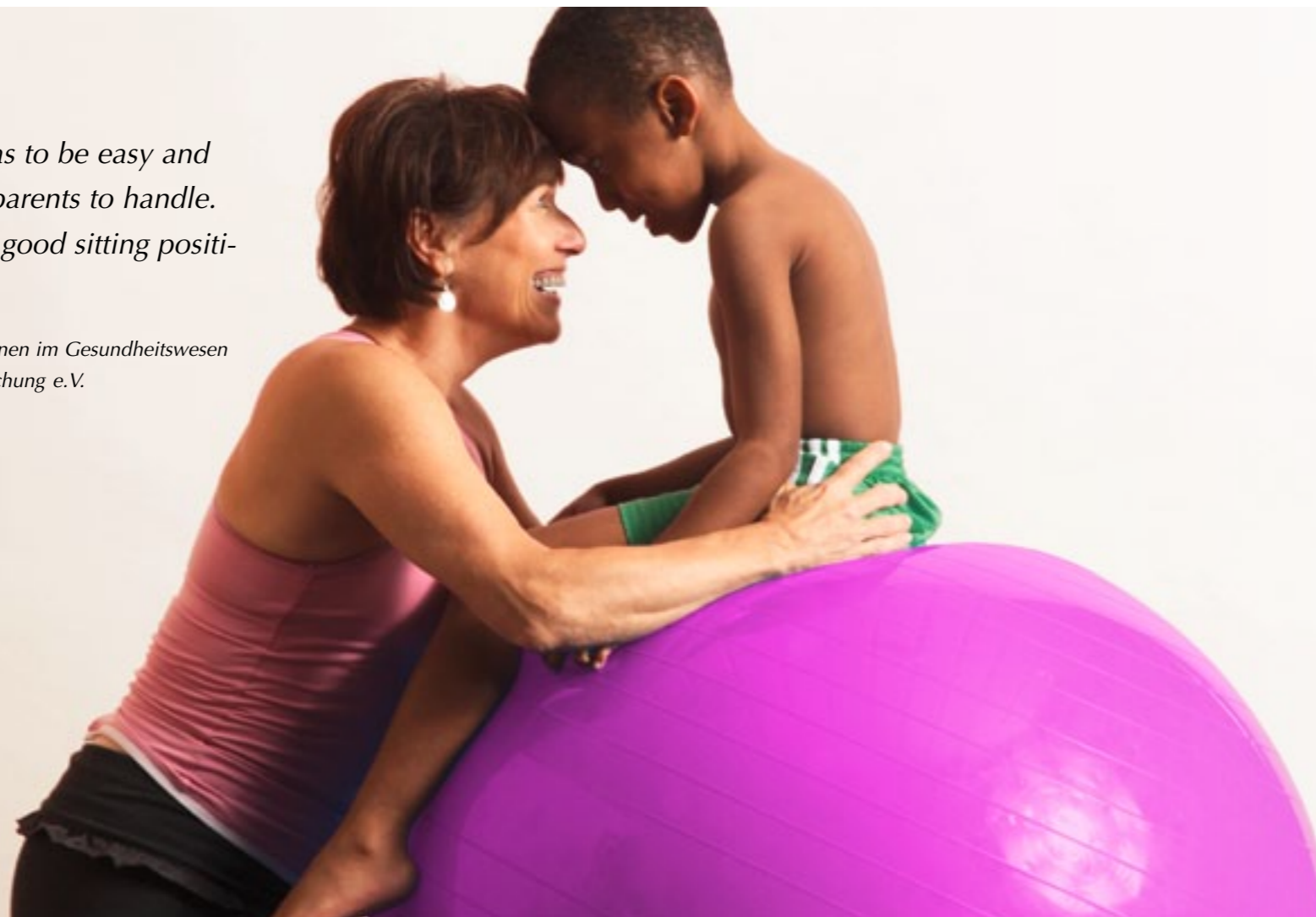


Easy, hassle-free and therapeutic adjustment for all involved

- Therapeutic adjustment options for an ideal sitting position
- Wide range of accessories for need-based positioning
- More child and less rehab stroller mean ideal integration

tRide – Rehab Stroller



The three-dimensionally adjustable headrest allows the child to sit relaxed and without tiring. It promotes an anatomically correct head and sitting position.

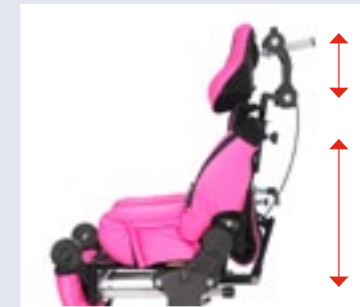

The short back and the separate head rest provide more freedom of movement from the torso to the upper extremities. The child sits in a stable position and at the same time can take part in his/her environment better.



The functional footrest adjusts to meet the child's individual needs. The stroller is easy to get into and out of which is good for muscle tone and promotes a posture that is easy on the back.



The stroller's sophisticated ergonomics allow it to adapt as the child grows offering perfect support over an extended period of time.



The sitting/reclined position is easily and quickly adjusted to the ideal position. tRide offers immediate and uncomplicated relief for the child's spine.



Technical Data

tRide Seat 1

tRide Seat 2

Seat depth	17 - 28 cm / 6.7 - 11"	24 - 35 cm / 9.4 - 13.8"
Seat width	20 - 30 cm / 7.9 - 11.8"	26 - 35 cm / 10.2 - 13.8"
Seat height *	26 - 42 cm / 10.2 - 16.5"	32 - 52 cm / 12.6 - 20.5"
Lower leg length	17 - 30 cm / 6.7 - 11.8"	25 - 36 cm / 9.8 - 14.2"
Hip angle	90° to 140°	90° to 140°
Knee angle	90° to 180°	90° to 180°
Foot angle	78° to 100°	78° to 100°
Height adjustment of arm support	10 - 14 cm / 3.9 - 5.5"	12 - 18 cm / 4.7 - 7.1"
Size of arm support (W x H)	21 x 10 cm / 8.3 x 3.9"	29 x 12 cm / 11.4 x 4.7"
Size of footrest (W x D)	29 x 18 cm / 11.4 x 7.1"	33 x 21 cm / 13 x 8.3"
Overall dimensions (L x W x H)	85 x 41 x 57 cm / 33.5 x 16.1 x 22.4"	98 x 46 x 46 cm / 38.6 x 18.1 x 18.1"
Folded dimensions without footrest	63 x 36 x 38 cm / 24.8 x 14.2 x 15"	76 x 46 x 50 cm / 30 x 18.1 x 19.7"
Weight of the seat *	8 kg / 17.6 lbs	10, 8 kg / 23.8 lbs
Max. carrying capacity	35 kg / 77.2 lbs	40 kg / 88.2 lbs

* without headrest

Technical Data

tRide Chassis 1

tRide Chassis 2

Seat tilt-in-space	-10° to +35°	-10° to +35°
Seat height	46 cm / 18.1"	47 cm / 18.5"
Push bar height	76 - 117 cm / 29.9 - 46.1"	81 - 120 cm / 31.9 - 47.2"
Overall dimensions (L x W x H)	75 x 60 x 77 cm / 29.5 x 23.6 x 30.3"	81 x 73 x 79 cm / 31.9 x 28.7 x 31.1"
Folded dimensions (L x W x H)	79 x 60 x 38 cm / 31.1 x 23.6 x 15"	83 x 67 x 38 cm / 32.7 x 26.4 x 15"
Swivel wheels size (front)	7,5"	7,5"
Wheel size (back)	10"	10"
Weight (without seat unit)	8,2 kg / 18.1 lbs	9 kg / 19.8 lbs
Max. carrying capacity (incl. seat)	45 kg / 99.2 lbs	50 kg / 110.2 lbs
Overall dimensions of the tRide rehabilitation stroller, ready for use	90 x 60 x 101 cm	110 x 67 x 99 cm
Folded dimensions of the tRide rehabilitation Stroller (L x B x H)	35.4 x 23.6 x 39.8"	43.3 x 26.4 x 39"
	81 x 60 x 63 cm	79 x 64 x 72 cm
	31.9 x 23.6 x 24.8"	31.1 x 25.2 x 28.3